

Humanization Design

# Digital display speed control knob


9 levels speed regulation, better control of the pumping speed,  
improve the success rate of the experiment



Easily Data Processing

# Running in low noise

Running noise  $\leq 40\text{DB}$ , low noise and no pollution

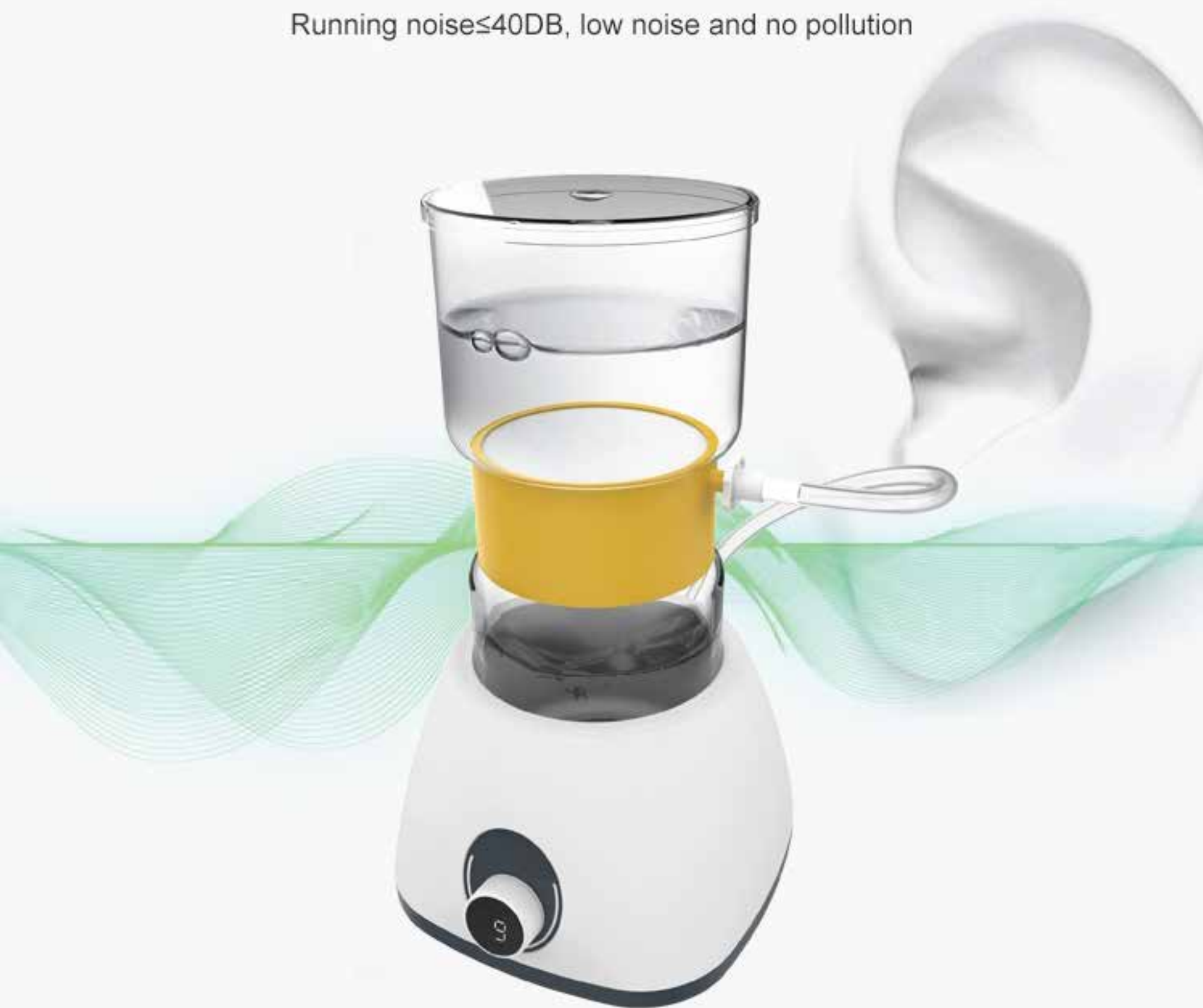

More Effective

# Designed for low power pumping

With adjustable pumping rate, it can be used with a disposable plastic vacuum filtration system.

Suitable for working scenarios with small processing capacity and slow pumping rate.

(Note: Large vacuum pump power is large, pumping rate is fixed, which can easily damage the plastic vacuum filtration system)



Not easy to damage  
the plastic device



Large capacity device can  
be connected externally

# ACCESSORIES



Machine



Power cord



Filter membrane



Silicone tube



使用说明书

User manual

## 6L Mini Vacuum Pump



9-speed speed regulation



Digital display speed control knob



Mini portable, strong power



Steady motion

### TECHNICAL DATA

Model	VP-6L	Pumping Rate	6L / Min
Type	Negative pressure suction	Ultimate Vacuum	0.06Mpa
Control Type	Digital display knob · 9 speeds	Power	10W